


*Guard columns  
WITHOUT press-  
tight connections-*

*protecting your  
analytical column has  
never been this easy!*

**•1 Innovative Integra-Guard Columns**

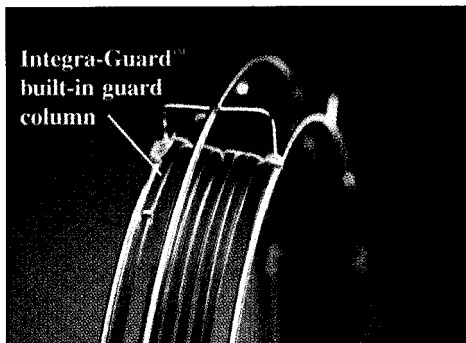

How many times has it been said about press-fit connectors, "Some people swear by them.. others swear at them"? For many analysts, the art of attaching a guard column to an analytical column is a mystery. Restek's chemists have discovered the solution to this mystery-the most reliable connection is no connection at all! No guard column system is more permanent than one continuous length of tubing containing both the guard column and the analytical column.

Restek now offers a wide variety of Integra-Guard™ capillary columns with a guaranteed leak-free connection! The column is tied using high temperature string with the guard column set apart from the analytical column. The transition area between the columns is the point at which the guard column ends and the analytical column begins. The entire setup is suspended in our unique "crush-free" cage that prevents the column from coming in contact with anything that could damage it.

The Integra-Guard column is so economical that we challenge you to beat our price with that of a conventional connection, even if you assemble it yourself. If you are currently using a guard column or considering one for the future, call your local distributor today and ask about our Integra-Guard columns.

*Ordering is easy. Just add the appropriate suffix number and price to the analytical column's catalog number and price. For example, a 30m, 0.25mm ID, 0.25um Rt@-5 with a 5-meter Integra-Guard column is cat.# 10223-124.*

mm ID	Length	Suffer #
0.25	5m	-124
	10m	-127
0.32	5m	-125
	10m	-128
0.53	5m	-126
	10m	-129

***Restek-bringing you  
breakthroughs in  
column protection!***

**Don McCandless,**  
Marketing Director and  
5th degree black belt in  
Isshinryu karate.

